

BLADDER MANAGEMENT

Liquick X-treme and Liquick X-treme Plus

Pre-activated, ready-to-use intermittent catheter systems



Teleflex

Liquick X-treme catheter system

Enjoy your independence – fast and simple

At Teleflex, we understand how important it is to help people make their lives as simple and as comfortable as possible, while assisting people to keep their independence and dignity. With this in mind, we have developed an innovative intermittent catheter system: the Liquick X-treme catheter system is a ready-to-use catheter for intermittent catheterisation that is designed to be intuitive and simple to apply virtually anywhere you need – for X-treme independence and comfort.





X-treme Speed

Ready-to-use

The Liquick X-treme catheter system is ready-to-use straight out of the packaging. There is no need to burst a sachet, to add lubrication or to wait for the hydrophilic coating to be activated. Due to the pre-activated hydrophilic coated SafetyCat catheter, the use of the Liquick X-treme catheter system is fast and easy.

X-treme Comfort

Simple non-touch handling

The modern packaging in an elegant design is easy to open. Equipped with a blue introduction aid and a sterile protective inner sleeve, which covers the complete surface of the catheter, the Liquick X-treme catheter system is designed to allow a comfortable non-touch catheterisation: all that you have to do is, open, remove, insert.

X-treme Convenience


Centrepiece of our technology

The Liquick X-treme catheter system is equipped with our SafetyCat catheter featuring softly rounded Soft Cat eyes, a flexible Ergothan tip as well as a pre-activated hydrophilic coating, the Liquick X-treme catheter system provides a gentle and smooth catheterisation.

Liquick X-treme Plus catheter system

All in one - suitable for daily use

The new ready-to-use Liquick X-treme Plus catheter system is a complete system for a quick, convenient catheterisation virtually at any place and any time. This catheter system has a preconnected urine collection bag making it convenient to use whenever you need to catheterise yourself. Additionally, all the features of the Liquick X-treme catheter system make this product suitable for daily use.

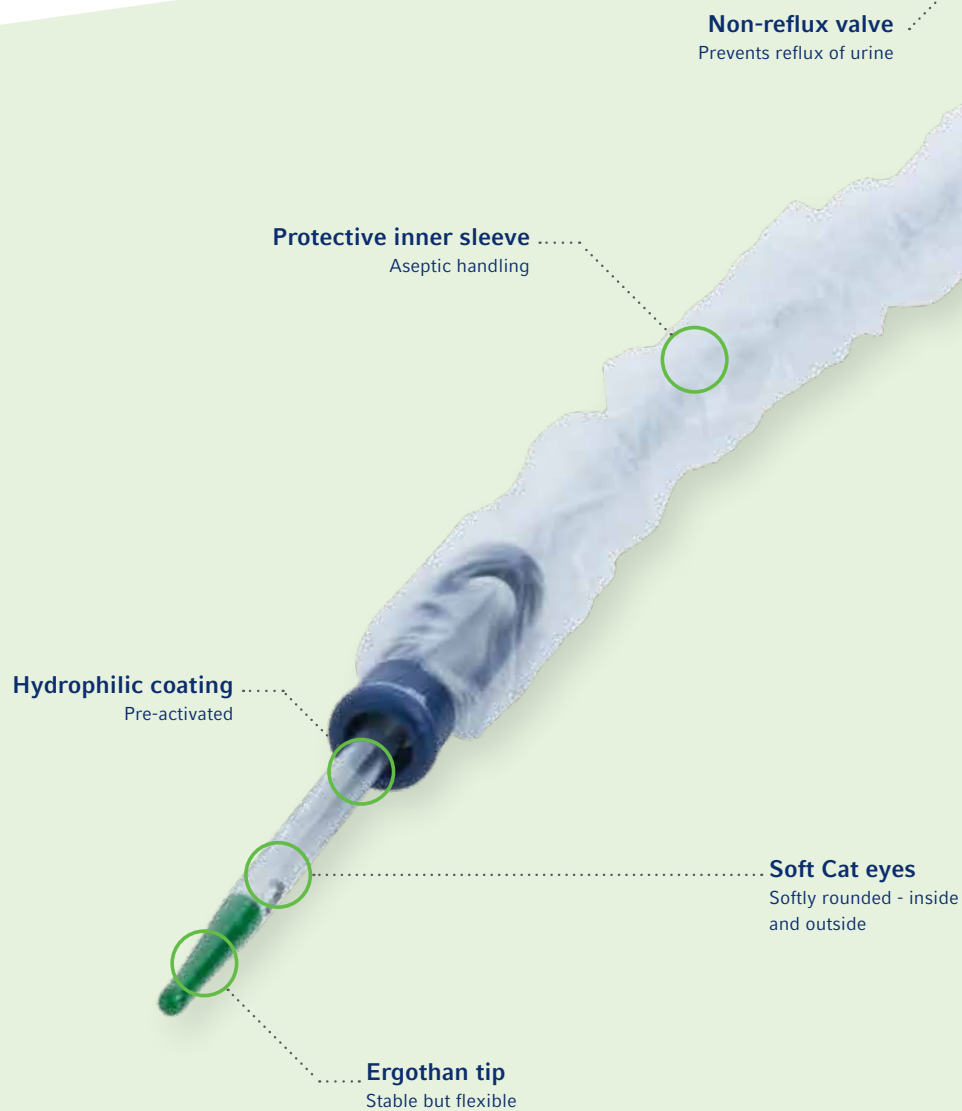



Perforation line
Easy emptying

Thumb hole
Facilitates handling

X-treme independence

The Liquick X-treme Plus catheter system is ready-to-use straight out of the packaging due to the pre-activated hydrophilic coated SafetyCat catheter. Equipped with a blue introduction aid, a protective inner sleeve and a preconnected urine collection bag this system is designed to allow an independent completely aseptic, non-touch catheterisation. Just open the packaging, remove the system, unfold the urine collection bag and insert the catheter.

Non-touch handling

Equipped with a blue introduction aid and a protective inner sleeve, which covers the complete surface of the catheter.

SafetyCat catheter

The Liquick X-treme Plus catheter system is equipped with the SafetyCat catheter with softly rounded Soft Cat eyes, a flexible Ergothan tip as well as a pre-activated hydrophilic coating designed for a gentle and smooth catheterisation.

Liquick X-treme catheter system



ERGOTHAN TIP – 2 EYES, STAGGERED

LENGTH	CH	REF	QTY
20 cm	10 ●	851822-000100	30
	12 ○	851822-000120	
	14 ●	851822-000140	
40 cm	10 ●	851842-000100	30
	12 ○	851842-000120	
	14 ●	851842-000140	
	16 ●	851842-000160	

Liquick X-treme Plus catheter system



ERGOTHAN TIP – 2 EYES, STAGGERED

LENGTH	CH	REF	QTY
20 cm	12 ○	852822-000120	30
	14 ●	852822-000140	
40 cm	12 ○	852842-000120	30
	14 ●	852842-000140	
	16 ●	852842-000160	



If you would like, you can order samples of the Liquick X-treme catheter system today.
Just scan the QR-code to order samples through our web site
www.teleflexurology.com



SafetyCat catheter

A solid base

All our catheter systems are equipped with the SafetyCat catheter.

Ergothan tip

The flexible Ergothan tip is designed to smoothly pass through the urethra and glide easily and comfortably into the bladder.



Soft Cat eyes

The Soft Cat eyes (drainage openings) are polished both on the inside and outside of the catheter, designed to allow gentle passage past the sensitive urethral mucosa.



Teleflex Headquarters International, Ireland

Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland
Phone +353 (0)9 06 46 08 00 · Fax +353 (0)14 37 07 73 · orders.intl@teleflex.com · www.teleflex.com

Teleflex

Grosvenor House, Horseshoe Crescent, Old Beaconsfield, Buckinghamshire, HP9 1LJ
Phone: +44 (0)1494 53 27 61 · Fax: +44 (0)1494 52 46 50

www.teleflexurology.com

Teleflex, the Teleflex logo, Liquick and Liquick X-treme are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Other names may be trademarks of their respective owners. Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised 03/21.

© 2022 Teleflex Incorporated. All rights reserved.

MCI-2021-0033-EN · REV 0 · 02 22 PDF

Please scan to visit
our website

**Teleflex**