



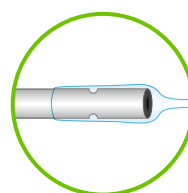

Rüsch  
**Superglide**  
Valve Ureteral Stent

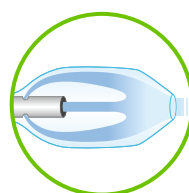

# Rüsch Superglide Valve Ureteral Stents

A stent designed to reduce vesicorenal and/or vesicoureteral reflux<sup>1,3</sup> and help enhance patient comfort over long term use compared to ureteral stents without anti-reflux valves.<sup>1,2,4</sup>




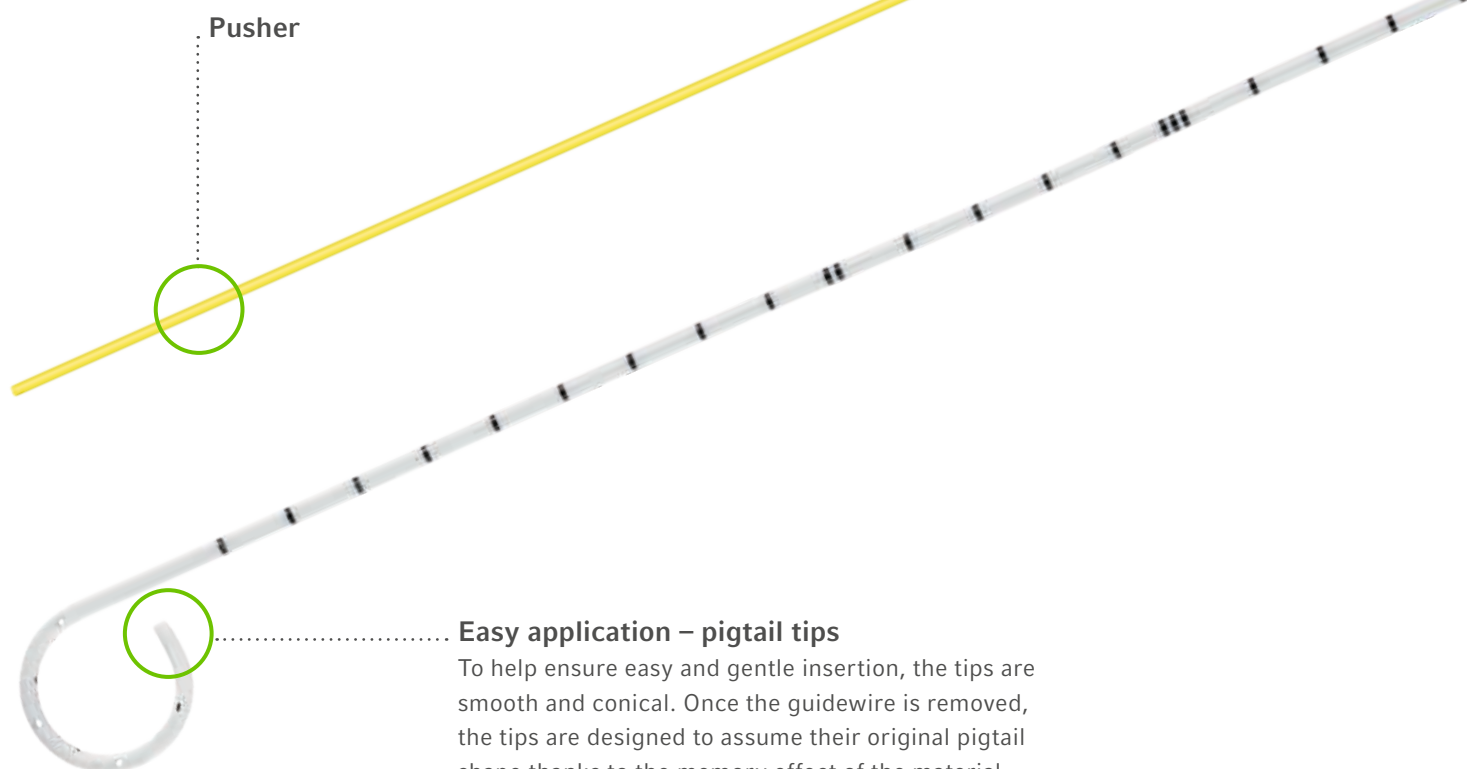
## Superglide Valve Ureteral Stents

Made of polyurethane with hydrogel coating

- Graduation in cm
- Stent completely radiopaque
- Pigtail and drainage eyes as indicated in the table
- Sterile, single use
- Not made with natural rubber latex

## Superglide Valve Ureteral Stent


TIP CONFIGURATION	EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART			ORDER SIZE (CH.)	QTY
		24 CM	26 CM	28 CM		
Central opening		334324	334326	334328	4.8/6/7/8	5



Pusher


### Easy application – pigtail tips

To help ensure easy and gentle insertion, the tips are smooth and conical. Once the guidewire is removed, the tips are designed to assume their original pigtail shape thanks to the memory effect of the material, which helps ensure secure retention in the renal pelvis.

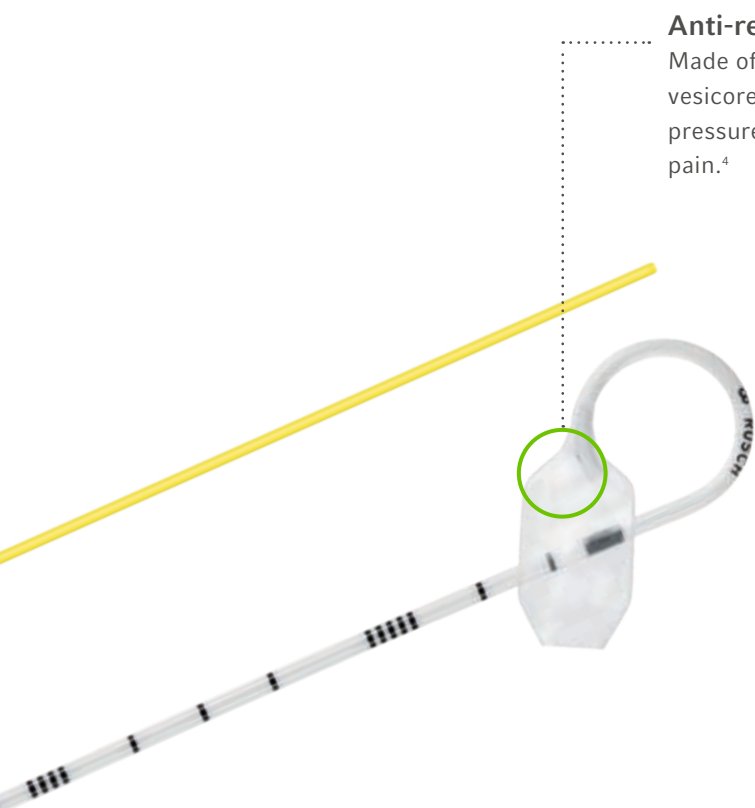

### Anti-reflux valve

Made of polyurethane at the bladder end of the stent designed to reduce vesicorenal and/or vesicoureteral reflux.<sup>3</sup> The valve is designed to avoid pressure increase in the upper urinary tract, which may help reduce flank pain.<sup>4</sup>

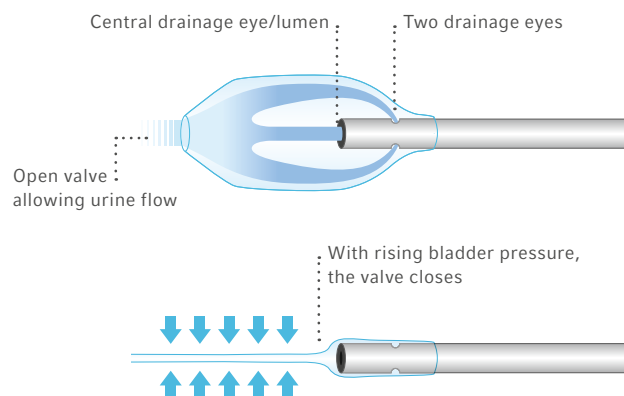

### Superglide Valve Ureteral Stent

**RÜSCH**

REF	EXCES CORRUM	QTY
334324-000048	Superglide anti-reflux valve ureteral stent, 4.8Ch, 24cm	5
334324-000060	Superglide anti-reflux valve ureteral stent, 6Ch, 24cm	5
334324-000070	Superglide anti-reflux valve ureteral stent, 7Ch, 24cm	5
334324-000080	Superglide anti-reflux valve ureteral stent, 8Ch, 24cm	5
334326-000048	Superglide anti-reflux valve ureteral stent, 4.8Ch, 26cm	5
334326-000060	Superglide anti-reflux valve ureteral stent, 6Ch, 26cm	5
334326-000070	Superglide anti-reflux valve ureteral stent, 7Ch, 26cm	5
334326-000080	Superglide anti-reflux valve ureteral stent, 8Ch, 26cm	5
334328-000048	Superglide anti-reflux valve ureteral stent, 4.8Ch, 28cm	5
334328-000060	Superglide anti-reflux valve ureteral stent, 6Ch, 28cm	5
334328-000070	Superglide anti-reflux valve ureteral stent, 7Ch, 28cm	5
334328-000080	Superglide anti-reflux valve ureteral stent, 8Ch, 28cm	5

All stent sizes accept a 0.035" guidewire.  
For retrograde application only

### References:

1. Nakada S, Patel S, Richie JP, Collins KA. *Placement and management of indwelling ureteral stents*, 2018, EBSCO Health, Ipswich, MA, DynaMed
2. Data on file at Teleflex. Test data may not be indicative of clinical performance.
3. Friedrich U, Vetter R, Jorgensen TM, Muhr M. *Experimental and Clinical Experience with New Pyeloureteral Stents in Pediatric Urology*. 1997;11(6):431-439
4. Ecke T, Bartel P, Hallman S, Ruttloff J. *Evaluation of Symptoms and Patient's Comfort for JJ-Ureteral Stents With and Without Antireflux Membrane Valve*. Urology 2010;75:212-216.

Teleflex is a global provider of medical technologies designed to improve the health and quality of people's lives. We apply purpose-driven innovation – a relentless pursuit of identifying unmet clinical needs – to benefit patients and healthcare providers. Our portfolio is diverse, with solutions in the fields of vascular and interventional access, surgical, anesthesia, cardiac care, urology, emergency medicine and respiratory care. Teleflex employees worldwide are united in the understanding that what we do every day makes a difference. For more information, please visit [teleflex.com](http://teleflex.com).

Teleflex is the home of Arrow, Deknatel, Hudson RCI, LMA, Pilling, Rüsch, UroLift, and Weck – trusted brands united by a common sense of purpose.

### **Corporate Office**

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

### **Regional Offices**

**United States:** Phone +1 919 544 8000, Toll Free 866 246 6990, [cs@teleflex.com](mailto:cs@teleflex.com), 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

**Latin America:** Phone +1 919 433 4999, [la.cs@teleflex.com](mailto:la.cs@teleflex.com), 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

**International:** Phone +353 (0)9 06 46 08 00, [orders.intl@teleflex.com](mailto:orders.intl@teleflex.com), Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

**Australia** 1300 360 226

**Austria** +43 (0)1 402 47 72

**Belgium** +32 (0)2 333 24 60

**Canada** +1 (0)800 387 9699

**China** (Shanghai) +86 (0)21 6060 7000

**China** (Beijing) +86 (0)10 6418 5699

**Czech Republic** +42 (0)495 759 111

**France** +33 (0)5 62 18 79 40

**Germany** +49 (0)7151 691 3004

**Greece** +30 210 67 77 717

**India** +91 (0)80 4093 4790

**Italy** +39 0362 58911

**Japan** +81 (0)3 6632 3600

**Korea** +82 2 536 7550

**Mexico** +52 55 5002 3500

**Netherlands** +31 (0)88 00 215 00

**New Zealand** 0800 601 100

**Poland** +48 22 4624032

**Portugal** +351 225 491 051

**Singapore** (SEA non-direct sales countries) +65 6439 3000

**Slovak Republic** +421 (0)3377 254 28

**South Africa** +27 (0)11 807 4887

**Spain** +34 91 198 84 31

**Switzerland** +41 (0)31 818 40 90

**United Kingdom** +44 (0)1494 53 27 61

**For more information, please visit [teleflex.com](http://teleflex.com)**

Teleflex, the Teleflex logo, Arrow, Deknatel, Hudson RCI, LMA, Pilling, Rüsch, Superglide, UroLift, and Weck are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised: 11/2019.