

# Urology catalogue

For almost every application



# Urology

Colour code and tips – Balloon catheters. . . . .3

1

Balloon catheters . . . . .4

1.1

Standard 2-way drainage . . . . .4

1.2

Postoperative drainage, 2- and 3-way . . . . .7

1.3

Balloon catheters for children . . . . .13

1.4

Special catheters . . . . .14

1.5

Accessories . . . . .15

2

Bladder catheters without balloon . . . . .16

2.1

Catheters made of pvc . . . . .16

3

Urine collectors . . . . .18

3.1

Preconnected . . . . .18

3.2

Closed urinary drainage systems . . . . .19

3.3

Accessories . . . . .23

4

Ureteral catheters . . . . .24

4.1

Hydrophilic coated ureteral catheters for standard applications . . . . .24

4.2

Ureteral catheters for standard applications . . . . .25

4.3

Ureteral catheter for injection of contrast medium . . . . .27

4.4

Accessories . . . . .28

5

Ureteral stents. . . . .29

5.1

Hydrophilic coated ureteral stents . . . . .29

5.2

Ureteral stents . . . . .33

5.3

Rigid ureteral stents . . . . .36

5.4

Ureteral stents for children . . . . .37

5.5

Special ureteral stents . . . . .39

5.6

Eco ureteral stents . . . . .41

6

Nephrostomy . . . . .42

6.1

Puncture sets . . . . .42

6.2

Exchange sets . . . . .43

6.3

Catheters for percutaneous nephrostomy . . . . .45

7

Guidewires . . . . .46

8

Stone extractors . . . . .47

8.1

Flat wire baskets . . . . .47

8.2

Round wire baskets . . . . .48


9

Special products . . . . .49

9.1







Mechanical ureter dilators . . . . .49





Index of reference numbers. . . . .50

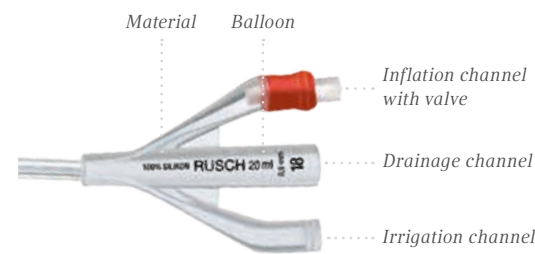


Information

## Teleflex balloon catheter tip variants and drainage eye configurations

TIPS	
Cylindrical	
Tiemann	
Couvellaire	
Dufour	
Mercier	
Whistle tip	

EYE ARRANGEMENTS			
			
1 eye	2 opposing eyes	2 eyes, side by side	2 eyes, staggered



Catheter external diameter, e.g. 18 Charrière/French = 6 mm

## Colour code

Ch. 6

Ch. 8

Ch. 10

Ch. 12

Ch. 14

Ch. 16

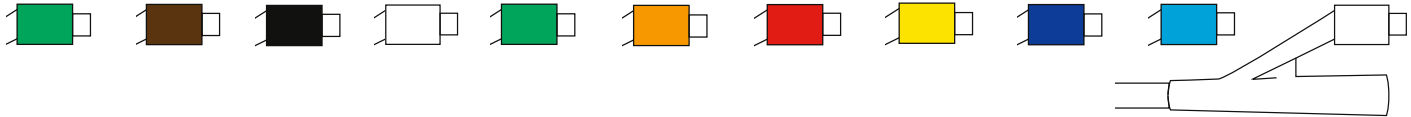
Ch. 18


Ch. 20

Ch. 22

Ch. 24

Ch. 26/28/30





**Contact us**  
Your requirements and suggestions are important to us. It is our ambition to find a solution that really meets your needs. You are welcome to contact us.

# 1 Balloon catheters

## 1.1 Standard 2-way drainage

### ProfilCath AquaFlate Glycerine Brilliant balloon catheter

Made of silicone, transparent, with longitudinal grooves designed for drainage of urethral secretion, with prefilled syringe for balloon inflation

- Prefilled syringe in sterile packaging (10 ml sterile 10 % glycerine solution)
- Valve for Luer and Luer-Lock syringe tips
- Radiopaque tip and contrast stripe along the catheter shaft
- 2 self-adhesive documentation labels
- Sterile, single use
- **Not made with natural rubber latex**



ProfilCath AquaFlate Glycerine						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	171505	<b>Cylindrical, solid tip</b>				
		Ch. 12–24	10 ml	2 opposing eyes	Approx. 41 cm	5



AquaFlate Glycerine						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	178000	<b>Cylindrical, solid tip</b>				
		Ch. 12–24	10 ml	2 opposing eyes	Approx. 40 cm	5
	178001	<b>For WOMEN, cylindrical, solid tip</b>				
		Ch. 12–24	10 ml	2 opposing eyes	Approx. 23 cm	5
	178305	<b>Tiemann</b>				
		Ch. 12–24	10 ml	1 eye	Approx. 41 cm	5



### Rüsch Brillant balloon catheter

Made of silicone, transparent

- Valve for Luer and Luer-Lock syringe tips
- Radiopaque tip and contrast stripe along the catheter shaft
- 2 self-adhesive documentation labels
- Sterile, single use
- **Not made with natural rubber latex**

Rüsch Brillant						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	170605	<b>Cylindrical, solid tip</b>				
		Ch. 12–24	10 ml	2 opposing eyes	Approx. 40 cm	10
	171305	<b>Tiemann</b>				
		Ch. 12–24	10 ml	1 eye	Approx. 41 cm	5

When ordering, please always indicate Ref. and Order size



Rüsch SiLaSil						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	186005*	<b>Cylindrical</b>				
		Ch. 12/14 Ch. 16–30	5–10 ml 5–15 ml	2 opposing eyes	Approx. 40 cm	10
	186505*	<b>Tiemann</b>				
		Ch. 12/14 Ch. 16–24	5–10 ml 5–15 ml	1 eye	Approx. 40 cm	10



Rüsch SympaCath						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	180505*	<b>Cylindrical</b>				
		Ch. 12/14 Ch. 16–30	5–10 ml 5–15 ml	2 opposing eyes	Approx. 40 cm	10



Rüsch Gold						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	180605*	<b>Cylindrical</b>				
		Ch. 12/14 Ch. 16–30	5–10 ml 5–15 ml	2 opposing eyes	Approx. 40 cm	10
	180105*	<b>For WOMEN, cylindrical</b>				
		Ch. 12/14 Ch. 16–22	5–10 ml 5–15 ml	2 opposing eyes	Approx. 25 cm	10
	181305*	<b>Tiemann</b>				
		Ch. 12/14 Ch. 16–26	5–10 ml 5–15 ml	1 eye	Approx. 40 cm	10

\* These products contain natural rubber latex which may cause allergic reactions.

### Rüsch SiLaSil balloon catheter

Made of latex, silicone elastomer inner and outer coating, blue

- Silicone elastomer coating on the inside and outside
- Valve for Luer and Luer-Lock syringe tips
- Sterile, single use

### Rüsch SympaCath balloon catheter

Made of latex, with hydrophilic coating, yellow

- Hydrophilic coating
- Valve for Luer and Luer-Lock syringe tips
- Sterile, single use

### Rüsch Gold balloon catheter

Made of latex, yellow

- Valve for Luer and Luer-Lock syringe tips
- For short-term drainage
- Siliconised
- Sterile, single use

SoftSimplastic balloon catheter

Made of PVC, transparent

- Latex balloon
- Valve for Luer and Luer-Lock syringe tips
- **Red** contrast stripe along the catheter shaft
- Sterile, single use



SoftSimplastic

	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	660610*	<b>Cylindrical</b>				
		Ch. 12–24	10 ml	2 opposing eyes	Approx. 42 cm	10
	661610*	<b>Tiemann</b>				
		Ch. 14–22	10 ml	1 eye	Approx. 42 cm	10

Rüsch balloon catheter

Made of polyurethane, transparent

- With Rüsch Gold latex balloon
- Valve for Luer and Luer-Lock syringe tips
- Sterile, single use



Wiruthan

	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	331305*	<b>Tiemann</b>				
		Ch. 6/8 Ch. 10/12	3 ml 5 ml	1 eye	Approx. 40 cm	5

1.2 Postoperative drainage, 2- and 3-way



Rüsch Brilliant

	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	170630	<b>Cylindrical, solid tip</b>				
		Ch. 18/20 Ch. 22/24	20 ml 30 ml	2 eyes, staggered	Approx. 40 cm	10
	172530	<b>Couvelaire</b>				
		Ch. 18/20 Ch. 22/24	20 ml 30 ml	2 eyes, staggered	Approx. 40 cm	5

Rüsch Brilliant balloon catheter

Made of silicone, transparent

- catheter for postoperative drainage, 2-way
- Valve for Luer and Luer-Lock syringe tips
  - Radiopaque tip and contrast stripe along the catheter shaft
  - 2 self-adhesive documentation labels
  - Sterile, single use
  - **Not made with natural rubber latex**

SoftSimplastic

	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	660630*	<b>Cylindrical</b>				
		Ch. 14–24	30 ml	2 opposing eyes	Approx. 42 cm	10
	661630*	<b>Tiemann</b>				
		Ch. 14–22	30 ml	1 eye	Approx. 42 cm	10
	662530* 662575*	<b>Couvelaire</b>				
		Ch. 16–24 Ch. 18–24	30 ml 75 ml	2 opposing eyes	Approx. 42 cm	10
	662630* 662675*	<b>Dufour</b>				
		Ch. 18–24 Ch. 18–24	30 ml 75 ml	2 eyes, staggered	Approx. 42 cm	10
	662430*	<b>Mercier</b>				
		Ch. 14–24	30 ml	2 eyes, staggered	Approx. 42 cm	10
	663630*	<b>Whistle tip</b>				
		Ch. 14–24	30 ml	1 eye	Approx. 42 cm	10

SoftSimplastic balloon catheter

Made of PVC, transparent

- Catheter for postoperative drainage, 2-way
- Latex balloon
  - Valve for Luer and Luer-Lock syringe tips
  - **Red** contrast stripe along the catheter shaft
  - Sterile, single use



When ordering, please always indicate Ref. and Order size

\* These products contain natural rubber latex which may cause allergic reactions.




Rüsch Gold balloon catheter

- Made of latex, yellow
- Catheter for postoperative drainage, 2-way
- Valve for Luer and Luer-Lock syringe tips
  - Siliconised
  - Sterile, single use

Rüsch Gold						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	180630*	Cylindrical				
		Ch. 12/14 Ch. 16–30	30 ml 30–50 ml	2 opposing eyes	Approx. 40 cm	10
	181330*	Tiemann				
		Ch. 14 Ch. 16–24	30 ml 30–50 ml	1 eye	Approx. 40 cm	10

Rüsch Gold balloon catheter





- Made of latex, armoured, yellow
- Catheter for postoperative drainage, 2-way
- Valve for Luer and Luer-Lock syringe tips
  - Siliconised
  - Sterile, single use

Rüsch Gold, Armoured						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	183530*	Cylindrical				
		Ch. 18–24	30–50 ml	2 eyes, staggered	Approx. 40 cm	10

The spiral in the wall enables the greatest flexibility possible, a consistently large lumen and high tension.

\* These products contain natural rubber latex which may cause allergic reactions.




Rüsch Brillant

	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	173830	Cylindrical, solid tip				
		Ch. 18/20 Ch. 22/24	20 ml 30 ml	2 eyes, staggered	Approx. 41 cm	5
	173810	Ch. 16/18 Ch. 20–24	80 ml 100 ml	2 eyes, staggered	Approx. 41 cm	5
	173930	Couvelaire				
		Ch. 18/20 Ch. 22/24	20 ml 30 ml	2 eyes, staggered	Approx. 41 cm	5
	173910	Ch. 16/18 Ch. 20–24	80 ml 100 ml	2 eyes, staggered	Approx. 41 cm	5
	174030	Dufour				
		Ch. 18/20 Ch. 22/24	20 ml 30 ml	2 eyes, staggered	Approx. 41 cm	5
	174130	Mercier				
		Ch. 18/20 Ch. 22/24	20 ml 30 ml	2 eyes, staggered	Approx. 41 cm	5

Rüsch Brillant balloon catheter

- Made of silicone, transparent
- Catheter for postoperative drainage, 3-way
- Valve for Luer and Luer-Lock syringe tips
  - Radiopaque tip and contrast stripe along the catheter shaft
  - 2 self-adhesive documentation labels (only Ref. 173830, 173930, 174030 and 174130)
  - Sterile, single use
  - **Not made with natural rubber latex**

ProfilCath

	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	173430	Cylindrical, solid tip				
		Ch. 18/20 Ch. 22/24	20 ml 30 ml	2 eyes, staggered	Approx. 41 cm	5
	173230	Couvelaire				
		Ch. 18/20 Ch. 22/24	20 ml 30 ml	2 eyes, staggered	Approx. 41 cm	5
	173330	With instillation lumen for medical treatment of the urethra, instillation opening below the balloon, cylindrical, solid tip				
		Ch. 18/20 Ch. 22/24	20 ml 30 ml	2 eyes, staggered	Approx. 41 cm	5

ProfilCath Rüsch Brillant balloon catheter

- Made of silicone, transparent, with longitudinal grooves designed for drainage of urethral secretion
- Catheter for postoperative drainage, 3-way
- Valve for Luer and Luer-Lock syringe tips
  - Radiopaque tip and contrast stripe along the catheter shaft
  - 2 self-adhesive documentation labels
  - Sterile, single use
  - **Not made with natural rubber latex**

When ordering, please always indicate Ref. and Order size

- Balloon catheter**  
Made of silicone, transparent  
Catheter for postoperative drainage, 3-way
- Large volume balloon fits closely to the shaft
  - Valve for Luer and Luer-Lock syringe tips
  - Radiopaque tip and contrast stripe along the catheter shaft
  - Sterile, single use
  - **Not made with natural rubber latex**



3-Way Silicone Foley						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	173800	<b>Cylindrical</b>				
		Ch. 16	40 ml	2 eyes, staggered	Approx. 42cm	5
		Ch. 18	50 ml			
		Ch. 20	60 ml			
		Ch. 22	70 ml			
		Ch. 24	80 ml			
	173900	<b>Couvelaire</b>				
		Ch. 16	40 ml	2 eyes, staggered	Approx. 42cm	5
		Ch. 18	50 ml			
		Ch. 20	60 ml			
		Ch. 22	70 ml			
		Ch. 24	80 ml			
	174000	<b>Dufour</b>				
		Ch. 16	40 ml	2 eyes, staggered	Approx. 42cm	5
		Ch. 18	50 ml			
		Ch. 20	60 ml			
		Ch. 22	70 ml			
		Ch. 24	80 ml			
	174100	<b>Mercier</b>				
		Ch. 16	40 ml	2 eyes, staggered	Approx. 42cm	5
		Ch. 18	50 ml			
		Ch. 20	60 ml			
		Ch. 22	70 ml			
		Ch. 24	80 ml			

SoftSimplastic


	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	664330*	<b>Cylindrical</b>				
		Ch. 18–24	30 ml	2 eyes, side by side	Approx. 42 cm	10
	664130* 664176*	<b>Couvelaire</b>				
		Ch. 16–24 Ch. 18–24	30 ml 75 ml	2 opposing eyes	Approx. 42 cm	10
	664530* 664575*	<b>Dufour</b>				
		Ch. 18–22 Ch. 18–24	30 ml 75 ml	2 eyes, staggered	Approx. 42 cm	10
	664230* 664275*	<b>Mercier</b>				
		Ch. 20–24 Ch. 18–24	30 ml 75 ml	2 eyes in front and 1 eye behind balloon 2 eyes, staggered	Approx. 42 cm	10
	663230*	<b>Whistle tip</b>				
		Ch. 18–24	30 ml	2 eyes, side by side	Approx. 42 cm	10






- SoftSimplastic balloon catheter**  
Made of PVC, transparent  
Catheter for postoperative drainage, 3-way
- Latex balloon
  - Valve for Luer and Luer-Lock syringe tips
  - **Red** contrast stripe along the catheter shaft
  - Sterile, single use

**Rüsch Gold balloon catheter**  
Made of latex, yellow  
Catheter for postoperative drainage, 3-way


- Valve for Luer and Luer-Lock syringe tips
- Siliconised
- Sterile, single use






Rüsch Gold						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	183405*	<b>Cylindrical</b>				
		Ch. 18–26	5–15 ml	2 eyes, side by side	Approx. 40 cm	10
	183430*	<b>Cylindrical</b>				
		Ch. 18–26	30–50 ml	2 eyes, side by side	Approx. 40 cm	10
	183230*	<b>Whistle tip</b>				
		Ch. 18–24	30–50 ml	2 eyes, side by side	Approx. 40 cm	10

**Rüsch Gold balloon catheter**  
Made of latex, armoured, yellow  
Catheter for postoperative drainage, 3-way

- Valve for Luer and Luer-Lock syringe tips
- Siliconised
- Sterile, single use






Rüsch Gold, Armoured						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	183830*	<b>Cylindrical</b>				
		Ch. 18–24	30–50 ml	2 eyes, side by side	Approx. 40 cm	10
	183960*	<b>Couvelaire</b>				
		Ch. 18–24	30–50 ml	2 eyes, staggered	Approx. 40 cm	10
	183930*	<b>Whistle tip</b>				
		Ch. 18–24	30–50 ml	2 eyes, side by side	Approx. 40 cm	10

The spiral in the wall enables the greatest flexibility possible, a consistently large lumen and high tension.

\* These products contain natural rubber latex which may cause allergic reactions.


### 1.3 Balloon catheters for children





AquaFlate Glycerine						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	178003	<b>Cylindrical, solid tip, with stylet</b>				
		Ch. 6 Ch. 8 Ch. 10	1.5 ml 3 ml 5 ml	2 opposing eyes	Approx. 31 cm	5
	178004	<b>Tiemann, with stylet</b>				
		Ch. 8 Ch. 10	3 ml 5 ml	1 eye	Approx. 31 cm	5

**AquaFlate Glycerine Brilliant balloon catheter**  
Made of silicone, transparent, with prefilled syringe for balloon inflation


- Prefilled syringe in sterile packaging (1.5, 3 or 5 ml sterile 10 % glycerine solution)
- Valve for Luer and Luer-Lock syringe tips
- Radiopaque tip and contrast stripe along the catheter shaft
- 2 self-adhesive documentation labels
- Sterile, single use
- **Not made with natural rubber latex**




Rüsch Brilliant						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	170003	<b>Cylindrical, solid tip, with stylet</b>				
		Ch. 6 Ch. 8 Ch. 10	1.5 ml 3 ml 5 ml	2 opposing eyes	Approx. 31 cm	5
	170004	<b>Tiemann, with stylet</b>				
		Ch. 8 Ch. 10	3 ml 5 ml	1 eye	Approx. 31 cm	5

**Rüsch Brilliant balloon catheter**  
Made of silicone, transparent

- Valve for Luer and Luer-Lock syringe tips
- Radiopaque tip and contrast stripe along the catheter shaft
- 2 self-adhesive documentation labels
- Sterile, single use
- **Not made with natural rubber latex**



Rüsch Gold						
	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	180003*	<b>Cylindrical, with stylet</b>				
		Ch. 8/10	3 ml	2 opposing eyes	Approx. 30 cm	10

**Rüsch Gold balloon catheter**  
Made of latex, yellow

- Valve for Luer and Luer-Lock syringe tips
- For short-term drainage
- Siliconised
- Sterile, single use

When ordering, please always indicate Ref. and Order size

## 1.4 Special catheters

### Balloon catheter for the neobladder for postoperative drainage

Made of silicone, transparent, with longitudinal grooves designed for drainage of secretion, 3-way

- Valve for Luer and Luer-Lock syringe tips
- Contrast stripe along the catheter shaft
- 2 self-adhesive documentation labels
- Sterile, single use
- **Not made with natural rubber latex**



#### Catheter for the neobladder

	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	173420	<b>Central open tip</b> , approx. 4 cm long				
		Ch. 20 Ch. 22	20 ml 30 ml	6 opposing eyes	Approx. 41 cm	5



### Supracath balloon catheter for suprapubic urine drainage

Made of silicone, transparent, 2-way

- Valve for Luer- and Luer-Lock syringe tips
- Radiopaque tip (only Ch. 14-20)
- Contrast stripe along the catheter shaft
- 2 self-adhesive documentation labels
- Sterile, single use
- **Not made with natural rubber latex**



#### SupraCath balloon catheter

	REF.	ORDER SIZE	BALLOON	EYES	LENGTH	QTY
	170606	<b>Central open tip</b>				
		Ch. 8 Ch. 10 Ch. 12–20	3 ml 5 ml 10 ml	2 opposing eyes	Approx. 41 cm	10

When ordering, please always indicate Ref. and Order size

## 1.5 Accessories

	REF.	DESCRIPTION	QTY
	556000	<b>Silkospray</b>	
		Silkospray is designed to avoid encrustation of rubber, latex and PVC devices.	10
	850550	<b>SilFlate</b> , prefilled syringe <sup>1</sup>	
		• With sterile 10 % aqueous glycerine solution • 10 ml • Sterile, single use	10

<sup>1</sup> Attention: To inflate silicone balloons only



## 2 Bladder catheters without balloon

### 2.1 Catheters made of PVC

**Red bladder catheter**

- Made of PVC
- Red
  - With funnel
  - Sterile, single use



Red bladder catheter					
	REF.	ORDER SIZE	EYES	LENGTH	QTY
	252000	<b>Mercier</b>			
		Ch. 8–24	1 eye	Approx. 38 cm	5
	259700	<b>Tiemann</b>			
		Ch. 6–24	1 eye	Approx. 38 cm	5

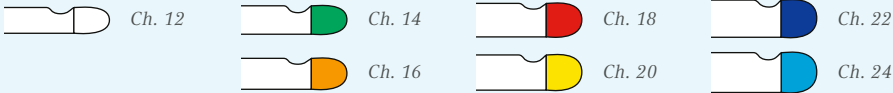
**Bladder catheter Tiemann**

- Made of PVC
- Transparent
  - With funnel
  - Sterile, single use



Bladder catheter Tiemann					
	REF.	ORDER SIZE	EYES	LENGTH	QTY
	221800	<b>Tiemann</b> , solid tip, siliconised, size acc. to colour code			
		Ch. 12–24	2 eyes, staggered	Approx. 40 cm	50

**Charrière Size is Colour Coded:**






When ordering, please always indicate Ref. and Order size



# 3 Urine collectors



## 3.1 Preconnected


Preconnected				
	REF.	COMPOSITION	ORDER SIZE	QTY
Sterile, single use, not made with natural rubber latex	171100 <sup>1</sup>	Preconnected		10
		<b>Rüsch Brilliant silicone balloon catheter</b> Cylindrical solid tip, 2 opposing eyes, length approx. 40 cm, balloon 10 ml	Ch. 12/14/16/18/20/22	
		<b>Rüsch S-Bag – closed urinary drainage system</b> Capacity 2,000 ml, length of tube 120 cm		
		<b>SilFlate – prefilled syringe<sup>2</sup></b> With sterile 10 % aqueous glycerine solution, 10 ml (Ref. 850550, if ordered separately)		
<b>Rüsch Brilliant silicone balloon catheter</b> <ul style="list-style-type: none"><li>• Transparent 2-way catheter with radiopaque tip and contrast stripe along the catheter shaft</li><li>• Cylindrical solid tip and 2 opposing eyes</li><li>• Enclosed prefilled syringe for inflation of the balloon, in sterile packaging, 10 ml sterile 10 % aqueous glycerine solution</li></ul> <b>Rüsch S-Bag</b> <ul style="list-style-type: none"><li>• Ridged universal connector</li><li>• Flat needle-free urine sampling port</li><li>• Vented<sup>3</sup> drip chamber with integrated anti-reflux valve</li><li>• Ventilation<sup>3</sup> at the bag</li><li>• Approximate graduation in ml</li><li>• Flexible inlet tubing with sheet clip and slide clamp</li><li>• Hook or cord options for fixation</li><li>• Drainage outlet with single hand clamp</li><li>• Doghouse holder to fix the outlet tube</li></ul>				

Connection between catheter and closed urinary drainage system

<sup>1</sup>Ref. 171100 comes with 2 self-adhesive documentation labels  
<sup>2</sup>Attention: To inflate silicone balloons only  
<sup>3</sup>Ventilation features a hydrophobic and oleophobic membrane


## 3.2 Closed urinary drainage systems

Rüsch S-Bag					
	REF.	NAME	CAPACITY	TUBE LENGTH	QTY
Sterile, single use	673044-000000 673049-000000 673049-002500	<b>Rüsch S-Bag</b> <b>Rüsch S-Bag A</b> <b>Rüsch S-Bag I</b>	2,000 ml 2,000 ml 2,500 ml	120 cm	20
	<ul style="list-style-type: none"><li>• Ridged universal connector with protective cap</li><li>• Flat needle-free urine sampling port</li><li>• Vented<sup>1</sup> drip chamber with integrated anti-reflux valve</li><li>• Ventilation<sup>1</sup> at the bag</li><li>• Approximate graduation in ml</li><li>• Flexible inlet tubing with sheet clip and slide clamp</li><li>• Hook or cord options for fixation</li><li>• Drainage outlet with single hand clamp (Ref. 673044-000000) or with T-tap valve (Ref. 673049-000000 and Ref. 673049-002500)</li><li>• Doghouse holder to fix the outlet tube</li></ul>			 T-tap outlet valve Ref. 673049-000000 and Ref. 673049-002500	


Rüsch P-Bag					
	REF.	NAME	CAPACITY	TUBE LENGTH	QTY
Sterile, single use	673052-000000	<b>Rüsch P-Bag</b>	1,000 ml	90 cm	20
	<ul style="list-style-type: none"><li>• Ridged universal connector with protective cap</li><li>• Flat needle-free urine sampling port</li><li>• Vented<sup>1</sup> drip chamber with integrated anti-reflux valve</li><li>• Ventilation<sup>1</sup> at the bag</li><li>• Approximate graduation in ml</li><li>• Flexible inlet tubing with sheet clip and slide clamp</li><li>• Hook or cord options for fixation</li><li>• Drainage outlet with single hand clamp</li><li>• Doghouse holder to fix the outlet tube</li></ul>				

<sup>1</sup> Ventilation features a hydrophobic and oleophobic membrane

Urine drainage systems consisting of urine meter and urine collector


Rüsch closed urine meter systems					
	REF.	CAPACITY URINE METER	CAPACITY URINE COLLECTOR	TUBE LENGTH	QTY
Sterile, single use	673081-000000 673082-000000	500 ml	2,000 ml	150 cm 120 cm	10
<div><div></div><div><p><b>Urine meter</b></p><ul style="list-style-type: none"><li>• 500 ml capacity with integrated measurement chamber (1 ml graduation)</li><li>• Ridged universal connector with protective cap</li><li>• Flat needle-free urine sampling port</li><li>• Tube reinforcement above drip chamber</li><li>• Vented<sup>1</sup> drip chamber with integrated anti-reflux valve</li><li>• Ventilation<sup>1</sup></li><li>• Flexible inlet tubing with sheet clip and slide clamp</li><li>• Fixation straps</li><li>• Drainage outlet</li><li>• Overflow inside the urine meter</li></ul><p><b>Urine collector</b></p><ul style="list-style-type: none"><li>• 2,000 ml capacity</li><li>• Anti-reflux valve</li><li>• Ventilation<sup>1</sup></li><li>• Approximate graduation in ml</li><li>• T-tap outlet valve</li><li>• Slit to fix the outlet tube</li></ul></div></div>					

Closed urinary drainage systems with integrated 400 ml measurement chamber

Rüsch U-Bag					
	REF.	NAME	CAPACITY	TUBE LENGTH	QTY
Sterile, single use	673064-000000 673065-000000	Rüsch U-Bag	2,000 ml	120 cm 150 cm	10
<div><div></div><div><ul style="list-style-type: none"><li>• Ridged universal connector with protective cap</li><li>• Flat needle-free urine sampling port</li><li>• Vented<sup>1</sup> drip chamber with integrated anti-reflux valve</li><li>• 400 ml mounted pivoting double measurement chamber with integrated 50 ml measurement chamber and sampling port</li><li>• Ventilation<sup>1</sup> at the measurement chamber and bag</li><li>• Approximate graduation in ml</li><li>• Flexible inlet tubing with sheet clip and slide clamp</li><li>• Hook and loop strap with D ring for fixation</li><li>• Drainage outlet with single hand clamp</li><li>• Doghouse holder to fix the outlet tube</li></ul></div></div>					

<sup>1</sup> Ventilation features a hydrophobic and oleophobic membrane


Exchange collector for Ref. 673081-000000 and Ref. 673082-000000

	REF.	CAPACITY	QTY
Sterile, single use	673083-000000	2,000 ml	50
<div><div></div><div><ul style="list-style-type: none"><li>• Anti-reflux valve</li><li>• Ventilation<sup>1</sup></li><li>• Approximate graduation in ml</li><li>• T-tap outlet valve</li><li>• Slit to fix the outlet tube</li></ul></div></div>			


<sup>1</sup> Ventilation features a hydrophobic and oleophobic membrane



Closed urinary drainage system for postoperative bladder irrigation and drainage

Rüsch TUR-Bag					
	REF.	NAME	CAPACITY	TUBE LENGTH	QTY
Sterile, single use	673053-000000	Rüsch TUR-Bag	4,000ml	120 cm	20
<div><div></div><div><ul style="list-style-type: none"><li>• Ridged universal connector with protective cap</li><li>• Needle-free urine sampling port</li><li>• Vented¹ drip chamber with integrated anti-reflux valve</li><li>• Ventilation¹ at the bag</li><li>• Approximate graduation in ml</li><li>• Flexible inlet tubing with sheet clip and slide clamp</li><li>• Hook or cord options for fixation</li><li>• T-tap outlet valve</li><li>• Doghouse holder to fix the outlet tube</li></ul></div></div>					


Closed urinary drainage system for postoperative bladder irrigation and drainage

Rüsch TUR-Bag with bulb					
	REF.	NAME	CAPACITY	TUBE LENGTH	QTY
Sterile, single use	673054-000000	Rüsch TUR-Bag	4,000ml	120 cm	10
<div><div></div><div><ul style="list-style-type: none"><li>• Ridged universal connector with protective cap</li><li>• Flat needle-free urine sampling port</li><li>• Vented¹ drip chamber with integrated anti-reflux valve</li><li>• Ventilation¹ at the bag</li><li>• Approximate graduation in ml</li><li>• Flexible inlet tubing with suction bulb, sheet clip and 2 slide clamps</li><li>• Hook or cord options for fixation</li><li>• T-tap outlet valve</li><li>• Doghouse holder to fix the outlet tube</li></ul></div></div>					

¹Ventilation features a hydrophobic and oleophobic membrane


### 3.3 Accessories

Urine collection bag to empty urine drainage systems and for the disposal of urine

Urine collection bag with ridged universal connector			
	REF.	CAPACITY	QTY
Non-sterile, single use	673045-000000	2,000ml	250
<div><div></div><div><ul style="list-style-type: none"><li>• 2,000ml bag with approximate ml graduation</li><li>• Ridged universal connector with protective cap¹ (suitable for urine drainage systems with <b>flexible</b> outlet tube)</li><li>• Slide clamp at the inlet tubing to close the bag</li></ul><p>¹The urine collection bag Ref. 673045-000000 ideally fits the closed urinary drainage systems: Preconnected (Ref. 171100), Rüsch S-Bag (Ref. 673044-000000), Rüsch P-Bag (Ref. 673052-000000), Rüsch U-Bag (Ref. 673064-000000 &amp; 673065-000000) as well as other urine drainage systems with flexible outlet tubes.</p></div></div>			

The urine collection bag can be connected directly to the outlet tube of the urinary drainage system. An alternative to using open urine collection methods.

Urine collection bags to empty urine drainage systems and for the disposal of urine

Urine collection bags with female connector			
	REF.	CAPACITY	QTY
Non-sterile, single use	673045-002000	2,000 ml	250
	673045-002500	2,500 ml	25
<div><div></div><div><ul style="list-style-type: none"><li>• 2,000 or 2,500ml bag with approximate ml graduation</li><li>• Female connector, suitable for urine drainage systems with <b>T-tap outlet valve</b>²</li><li>• Slide clamp at the inlet tubing to close the bag</li></ul><p>²The urine collection bags Ref. 673045-002000 and Ref. 673045-002500 ideally fit the closed urinary drainage systems: Rüsch S-Bag A (Ref. 673049-000000), Rüsch S-Bag I (Ref. 673049-002500), Rüsch closed urine meter systems (Ref. 673081-000000 and Ref. 673082-000000) as well as other urine drainage systems with T-tap outlet valve.</p></div></div>			



**Information**

When being emptied, the urine drainage system forms a firm and closed unit with the urine collection bag. After emptying the urine drainage system the urine collection bag can be closed using the slide clamp and then disposed of.









# 4 Ureteral catheters

## 4.1 Hydrophilic coated ureteral catheters for standard applications

- Superglide ureteral catheters**
- Made of PVC with hydrophilic coating
  - Radiopaque
  - Graduation in cm
  - 2 different markers to identify right/left side
  - Enclosed Luer-Lock adapter
  - With stylet
  - Sterile, single use



Superglide Ureteral catheters					
TIP CONFIGURATION	REF.	ORDER SIZE	LENGTH	EYES	QTY
With straight tip					
 Cylindrical, closed	264701	Ch. 3–7	Approx. 70 cm	1	10
 Cylindrical, closed	264702	Ch. 3–6	Approx. 70 cm	2	10
For coaxial applications					
 Central opening	264710	Ch. 3–8	Approx. 70 cm	–	10
 Central opening	264711	Ch. 3–8	Approx. 70 cm	1	10
With curved tip					
For coaxial applications					
 Curved, central opening	264810	Ch. 3–8	Approx. 70 cm	–	10
 Curved, central opening	264811	Ch. 4–7	Approx. 70 cm	1	10

Maximum diameter of guidewire for Superglide ureteral catheters with central opening					
REF. 264710 / 264711 / 264810 / 264811					
Ch. 3	Ch. 4	Ch. 5	Ch. 6	Ch. 7	Ch. 8
0.018"	0.028"	0.035"	0.038"	0.038"	0.038"
0.45 mm	0.71 mm	0.89 mm	0.96 mm	0.96 mm	0.96 mm

## 4.2 Ureteral catheters for standard applications



Ureteral catheters					
TIP CONFIGURATION	REF.	ORDER SIZE	LENGTH	EYES	QTY
With straight tip					
 Cylindrical, closed	263602	Ch. 3–7	Approx. 70 cm	1	10
 Cylindrical, closed	223602	Ch. 3–8	Approx. 70 cm	1	10
 Cylindrical, closed*	223600	Ch. 3–7	Approx. 70 cm	1	10
 Cylindrical, closed	223606	Ch. 3–8	Approx. 70 cm	2	10
 Olive tip	263702	Ch. 3/4	Approx. 70 cm	1	10
 Whistle tip	263802	Ch. 4–7	Approx. 70 cm	1	10
For coaxial applications					
 Central opening	264600	Ch. 4/4.8	Approx. 70 cm	–	10
 Central opening	223604	Ch. 3–8	Approx. 70 cm	2	10

\*Soft PVC




Maximum diameter of guidewire for ureteral catheters with central opening					
REF. 264600		REF. 223604			
Ch. 4	Ch. 4.8	Ch. 3	Ch. 4/5	Ch. 6	Ch. 7/8
0.032"	0.038"	0.015"	0.018"	0.032"	0.045"
0.80 mm	0.96 mm	0.38 mm	0.45 mm	0.80 mm	1.14 mm

- Ureteral catheters**
- Made of PVC
  - Radiopaque
  - Graduation in cm
  - 2 different markers to identify right/left side
  - Enclosed Luer-Lock adapter
  - With stylet
  - Sterile, single use

When ordering, please always indicate Ref. and Order size

- Ureteral catheters
- Made of PVC
  - Radiopaque
  - Graduation in cm
  - 2 different markers to identify right/left side
  - Enclosed Luer-Lock adapter
  - With stylet
  - Sterile, single use




Ureteral catheters					
TIP CONFIGURATION	REF.	ORDER SIZE	LENGTH	EYES	QTY
With curved tip					
 Curved olive tip	264001	Ch. 3–7	Approx. 70 cm	1 20 mm from tip	10
For coaxial applications					
 Curved, central opening	264100	Ch. 3–8	Approx. 70 cm	–	10
 Curved, central opening	264102	Ch. 5–7	Approx. 70 cm	1	10

Maximum diameter of guidewire for Ureteral catheters with central opening			
REF. 264100			REF. 264102
Ch. 3	Ch. 4	Ch. 5/6/7/8	Ch. 5/6/7
0.015"	0.018"	0.035"	0.035"
0.38 mm	0.45 mm	0.89 mm	0.89 mm

4.3 Ureteral catheter for injection of contrast medium





Ureteral catheter										
TIP CONFIGURATION	REF.	ORDER SIZE			LENGTH		EYES	QTY		
 Chevassu, whistle tip	267500	Ureteral catheter, Chevassu pattern, for routine injection of contrast medium into the ureter during retrograde investigation								
		Ch. 4–10			Approx. 70 cm		–		10	
		Diameter of tip/shaft								
		Tip	Ch. 4	Ch. 5	Ch. 6	Ch. 7	Ch. 8	Ch. 9	Ch. 10	
		Shaft	Ch. 3	Ch. 3	Ch. 4	Ch. 4	Ch. 5	Ch. 6	Ch. 7	

- Ureteral catheter
- Made of PVC
  - Radiopaque
  - Graduation in cm
  - 2 different markers to identify right/left side
  - Enclosed Luer-Lock adapter
  - With stylet
  - Sterile, single use

Measuring Range Ref. 267500



4.4 Accessories

	REF.	DESCRIPTION	QTY
	340702	<b>Luer-Lock adapter</b>	10
		<ul style="list-style-type: none"><li>• Blue</li><li>• Clamped</li><li>• For Ch. 3–6</li><li>• Sterile, single use</li></ul>	
	340703	<b>Luer-Lock adapter</b>	10
		<ul style="list-style-type: none"><li>• White</li><li>• Clamped</li><li>• For Ch. 6–10</li><li>• Sterile, single use</li></ul>	

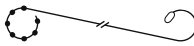
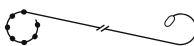


5 Ureteral stents

5.1 Hydrophilic coated ureteral stents



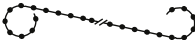
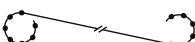
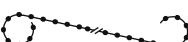
Superglide Valve DD ureter stents

TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART		ORDER SIZE (CH.)	QTY
	26 CM	28 CM		
Closed tip - for simultaneous applications				
	335340	335341	4.8/6/7	5
Closed tip - for coaxial applications				
	335440	335441	4.8/6/7/8	5

At the valve the O.D. of the stent is increased



Superglide DD ureter stents

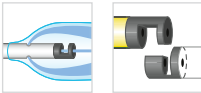
TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART				ORDER SIZE (CH.)	QTY
	22 CM	24 CM	26 CM	28 CM		
Closed tip - for simultaneous applications						
	–	334846	334840	–	4.8/6/7	5
	–	–	–	334841	4.8/6/7/8	
	–	–	334940	334941	4.8/6/7	
Open tip - for coaxial applications						
	335142	335143	–	–	4.8/6/7	5
	–	–	335140	335141	4.8/6/7/8	

Set Components

- Ureteral stent
- Fixation clamp
- Introducer made of polyurethane, yellow, approx. 45 cm long, with DD coupling
- Flexible guidewire made of stainless steel with PTFE coating and movable core: approx. 100 cm long in sets with cylindrical closed tip and approx. 150 cm long in sets with tip with central opening
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Maximum diameter of guidewire for Superglide Valve DD and Superglide DD ureteral stents with open tip

REF. 335440 / 335441 / 335142 / 335143 / 335140 / 335141			
Ch. 4.8	Ch. 6	Ch. 7	Ch. 8
0.028"	0.028"	0.035"	0.035"

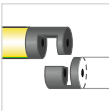


Superglide Valve DD ureteral stents

Made of polyurethane with hydrophilic coating, with DD coupling and anti-reflux valve at the vesical end

Directable with anti-reflux valve

- Graduation in cm
- Radiopaque stent
- Pigtail and drainage eyes as indicated in the table
- Supplied as a preassembled unit
- Sterile, single use



Superglide DD ureteral stents

Made of polyurethane with hydrophilic coating, with DD coupling at the vesical end

Directable

- Graduation in cm
- Radiopaque stent
- Pigtail and drainage eyes as indicated in the table
- Supplied as a preassembled unit
- Sterile, single use

When ordering, please always indicate Ref. and Order size


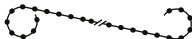
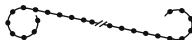
Superglide Integral ureteral stents

Made of polyurethane with hydrophilic coating, with retraction thread at the vesical end

- Graduation in cm
- Radiopaque stent
- Pigtail and drainage eyes as indicated in the table
- Sterile, single use



Superglide Integral ureter stents

TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART					ORDER SIZE (CH.)	QTY
	22 CM	24 CM	26 CM	28 CM	30 CM		
Closed tip – for simultaneous applications							
	–	334250	–	–	334249	4.8/6	5
	–	–	334247	334248	–	4.8/6/7	
Open tip – for coaxial applications							
	334742	–	–	–	–	6/7	5
	–	334743	–	–	–	4.8/6/7	
	–	–	334740	334741	–	4.8/6/7/8/9	
Open tip – for coaxial applications, without guidewire							
	–	–	334745	–	–	4.8/6/7	5
	–	334744	–	334746	–	6/7	


Set components

- Ureteral stent
- 2 fixation clamps
- Introducer made of polyurethane, yellow, approx. 45 cm long
- Flexible guidewire made of stainless steel with PTFE coating, fixed core, flexible safety tip and rigid tip: approx. 100 cm long in sets with cylindrical closed tip and approx. 150 cm long in sets with tip with central opening
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Maximum diameter of guidewire for Superglide Integral ureteral stents with open tip

REF. 334742 / 334743 / 334740 / 334741 / 334744 334745 / 334746				
Ch. 4.8	Ch. 6	Ch. 7	Ch. 8	Ch. 9
0.035"	0.035"	0.035"	0.038"	0.038"

Superglide Integral plus stents

TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART			ORDER SIZE (CH.)	QTY
	24 CM	26 CM	28 CM		
Open tip - for coaxial applications					
	–	334752	–	4.8/6/7	5
	–	–	334753	4.8/6/7/8	
	334751	–	–	6/7	

Set components


- Ureteral stent
- 2 fixation clamps
- Introducer made of polyurethane, yellow, approx. 45 cm long
- Flexible guidewire made of stainless steel with PTFE coating, fixed core, flexible safety tip and rigid tip: approx. 150 cm long
- Superglide ureteral catheter with central opening
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Maximum diameter of guidewire for Superglide Integral Plus ureteral stents with open tip

REF. 334751 / 334752 / 334753			
Ch. 4.8	Ch. 6	Ch. 7	Ch. 8
0.035"	0.035"	0.035"	0.038"



Superglide valve ureter stents

TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART			ORDER SIZE (CH.)	QTY
	24 CM	26 CM	28 CM		
Open tip – for coaxial applications, without guidewire					
	334324	334326	334328	4.8/6/7/8	5

At the valve the O.D. of the stent is increased.

Set components

- Ureteral stent
- Introducer made of polyurethane, yellow, approx. 45 cm long
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Superglide Integral ureteral stents

Made of polyurethane with hydrophilic coating, with retraction thread at the vesical end

With Superglide ureteral catheter, central opening, Ch. 5, max. I.D. 0.038"

- Graduation in cm
- Radiopaque stent
- Pigtail and drainage eyes as indicated in the table
- Sterile, single use



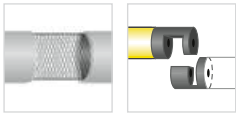
Superglide Valve ureter stents

Made of polyurethane with hydrophilic coating and valve at the vesical end

- Graduation in cm
- Radiopaque stent
- Pigtails and drainage eyes as indicated in the table
- Stents shown in this table accept a 0.035" guidewire
- Sterile, single use

When ordering, please always indicate Ref. and Order size





Superglide DD Tumour ureteral stents

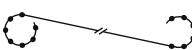
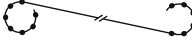
Made of polyurethane with hydrophilic coating, with DD coupling at the vesical end

Higher resistance to collapse when compared to the Teleflex standard polyurethane ureteral stents, due to a reinforced stent wall

Directable

- Graduation in cm
- Radiopaque stent
- Pigtail and drainage eyes as indicated in the table
- Supplied as a preassembled unit
- Sterile, single use



Superglide DD Tumour ureter stents							
TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART					ORDER SIZE (CH.)	QTY
	20 CM	22 CM	24 CM	26 CM	28 CM		
Closed tip – for simultaneous applications							
	–	–	335524	–	–	7	2
	–	–	–	335526	335528	6/7	
Open tip – for coaxial applications							
	335620	335622	335624	–	–	7	2
	–	–	–	335626	335628	6/7/8	

Set components

- Ureteral stent
- Fixation clamp
- Introducer made of polyurethane, yellow, approx. 45 cm long, with DD coupling
- Flexible guidewire made of stainless steel with PTFE coating and movable core: approx. 100 cm long in sets with cylindrical closed tip and approx. 150 cm long in sets with tip with central opening
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Maximum diameter of guidewire for Superglide DD Tumour ureteral stents with open tip		
REF. 335620 / 335622 / 335624 / 335626 / 335628		
Ch. 6	Ch. 7	Ch. 8
0.028"	0.035"	0.035"

5.2 Ureteral stents



Valve DD ureteral stents						
TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART				ORDER SIZE (CH.)	QTY
	22 CM	24 CM	26 CM	28 CM		
Closed tip - for simultaneous applications						
	335312 <sup>1</sup>	335311	335300	335301	4.8/6/7	5
Open tip - for coaxial applications						
	–	335411	–	–	4.8/6/7	5
	–	–	335400	335401	4.8/6/7/8	

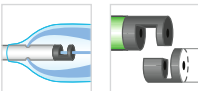
At the valve the O.D. of the stent is increased.

<sup>1</sup>Packaging unit of Ref. 335312 in Ch. 4.8 contains 2 stents

Set components

- Ureteral stent
- Fixation clamp
- Introducer made of polyurethane, green, approx. 45 cm long, with DD coupling
- Flexible guidewire made of stainless steel with PTFE coating and movable core: approx. 100 cm long in sets with cylindrical closed tip and approx. 150 cm long in sets with tip with central opening
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Maximum diameter of guidewire for Valve DD ureteral stents with open tip			
REF. 335411 / 335400 / 335401			
Ch. 4.8	Ch. 6	Ch. 7	Ch. 8
0.028"	0.028"	0.035"	0.035"

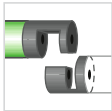


Valve DD ureteral stents

Made of polyurethane, with DD coupling and anti-reflux valve at the vesical end

Directable with anti-reflux valve

- Graduation in cm
- Radiopaque stent
- Pigtail and drainage eyes as indicated in the table
- Supplied as a preassembled unit
- Sterile, single use



- DD ureteral stents**
- Made of polyurethane, with DD coupling at the vesical end
- Directable
- Graduation in cm
  - Radiopaque stent
  - Pigtail and drainage eyes as indicated in the table
  - Supplied as a preassembled unit
  - Sterile, single use



**DD ureteral stents**

TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART					ORDER SIZE (CH.)	QTY
	22 CM	24 CM	26 CM	28 CM	30 CM		
Closed tip – for simultaneous applications							
	334812	–	–	–	–	4.8/6	5
	–	334811	–	–	334802	4.8/6/7	
	–	–	334800	–	–	4.8/6/7/8	
	–	–	–	334801	–	4.8/6/7/8/9	
	–	–	334900	334901	–	4.8/6/7	
Open tip – for coaxial applications							
	–	–	335100	–	335102	4.8/6/7/8	5
	335112	–	–	–	–	6/7	
	–	335111	–	–	–	6/7/8	
	–	–	–	335101	–	4.8/6/7/8/9	
	–	–	–	335166	–	4.8/6/7/8/9	
Closed tip – for ureterorenoscopes							
	–	–	335200	335201	–	4.8	5

- Set Components**
- Ureteral stent
  - Fixation clamp
  - Introducer made of polyurethane, green, approx. 45 cm long, with DD coupling (approx. 90 cm long in sets for ureterorenoscopes)
  - Flexible guidewire made of stainless steel with PTFE coating and movable core: approx. 100 cm long in sets with cylindrical closed tip and approx. 150 cm long in sets with tip with central opening and for ureterorenoscopes
  - Patient card
  - 3 self-adhesive documentation labels
  - Sterile packaging

**Maximum diameter of guidewire for DD ureteral stents with open tip**

335112 / 335111 / 335100 / 335101 / 335102 / 335166				
Ch. 4.8	Ch. 6	Ch. 7	Ch. 8	Ch. 9
0.028"	0.028"	0.035"	0.035"	0.035"



**Integral ureter stents**

TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART					ORDER SIZE (CH.)	QTY
	22 CM	24 CM	26 CM	28 CM	30 CM		
Closed tip – for simultaneous applications							
	334212	334211	–	–	–	4.8/6/7	5
	–	–	334200	334201	334202	4.8/6/7/8	
Open tip – for coaxial applications							
	334710	–	–	–	–	4.8/6/7	5
	–	334711	–	–	334714	4.8/6/7/8	
	–	–	334702	334703	–	4.8/6/7/8/9	
Closed tip – for ureterorenoscopes							
	–	–	334205	334206	–	4.8	5

- Integral ureteral stents**
- Made of polyurethane, with retraction thread at the vesical end
- Graduation in cm
  - Radiopaque stent
  - Pigtail and drainage eyes as indicated in the table
  - Sterile, single use

- Set Components**
- Ureteral stent
  - 2 fixation clamps
  - Introducer made of polyurethane, green, approx. 45 cm long (approx. 90 cm long in sets for ureterorenoscopes)
  - Flexible guidewire made of stainless steel with PTFE coating, fixed core, flexible safety tip and rigid tip: approx. 100 cm long in sets with cylindrical closed tip and approx. 150 cm long in sets with tip with central opening
  - For ureterorenoscopes: flexible guidewire made of stainless steel with PTFE coating and movable core, approx. 150 cm long
  - Patient card
  - 3 self-adhesive documentation labels
  - Sterile packaging

**Maximum diameter of guidewire for integral ureteral stents with open tip**

REF. 334710 / 334711 / 334702 / 334703 / 334714				
Ch. 4.8	Ch. 6	Ch. 7	Ch. 8	Ch. 9
0.035"	0.035"	0.035"	0.038"	0.038"

When ordering, please always indicate Ref. and Order size


5.3 Rigid ureteral stents

Integral ureteral stents

Made of firm polyurethane, with retraction thread at the vesical end

- Graduation in cm
- Radiopaque introducer and stent
- Pigtail and drainage eyes as indicated in the table
- Sterile, single use

Integral ureteral stents

TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART					ORDER SIZE (CH.)	QTY
	22 CM	24 CM	26 CM	28 CM	30 CM		
Open tip - for coaxial applications							
	334720	334721	334722	334723	334724	4.8/6/7/8/9	5

Set components

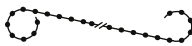
- Ureteral stent
- 2 fixation clamps
- Introducer made of firm polyurethane, green, approx. 45 cm long
- Flexible guidewire made of stainless steel with PTFE coating, with fixed core, flexible safety tip and rigid tip, approx. 150 cm long
- Patient card
- 3 documentation labels
- Sterile packaging

Superglide Integral ureteral stents

Made of firm polyurethane with hydrophilic coating and retraction thread at the vesical end

- Graduation in cm
- Radiopaque introducer and stent
- Pigtail and drainage eyes as indicated in the table
- Sterile, single use

Superglide Integral ureteral stents

TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART			ORDER SIZE (CH.)	QTY
	24 CM	26 CM	28 CM		
Open tip - for coaxial applications					
	334764	334765	334766	4.8/6/7/8	5

Set components

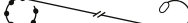
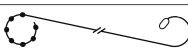
- Ureteral stent
- 2 fixation clamps
- Introducer made of firm polyurethane, yellow, approx. 45 cm long
- Patient card
- 3 documentation labels
- Sterile packaging

Maximum Diameter of Guidewire for Rigid Integral Stents with Central Opening

Ch. 4.8	Ch. 6	Ch. 7	Ch. 8	Ch. 9
0.035"	0.035"	0.035"	0.038"	0.038"

5.4 Ureteral stents for children

Valve DD ureteral stents

TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART						ORDER SIZE (CH.)	QTY
	10CM	12 CM	15 CM	18 CM	20 CM	22 CM		
Closed tip - for simultaneous applications								
	335315	335309	335314	335310	335313	335312	4.8	2
Open tip - for coaxial applications								
	–	335413	335414	335415	335416	335417	4.8	2

At the valve the O.D. of the stent is increased

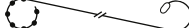
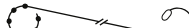
Set Components

- Ureteral stent
- Fixation clamp
- Introducer made of polyurethane, green, approx. 45 cm long, with DD coupling
- Flexible guidewire made of stainless steel with PTFE coating and movable core: approx. 100 cm long in sets with cylindrical closed tip and approx. 150 cm long in sets with tip with central opening
- Sets for coaxial application include an additional ureteral catheter with central opening in stent size (see table)
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Maximum diameter of guidewire for paediatric Valve DD Ureter Stents with open tip

REF. 335413 / 335414 / 335415 / 335416 / 335417	
Ch. 4.8	
0.028"	

Valve ureteral stents

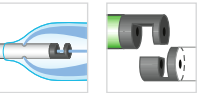
TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART						ORDER SIZE (CH.)	QTY
	10 CM	12 CM	15 CM	18 CM	20 CM	22 CM		
Closed tip - for simultaneous applications								
	–	335303	335304	335305	335306	335307	3/4	2
	335316	–	–	–	–	–	4	
Open tip - for coaxial applications								
	–	335403	335404	335405	335406	335407	3/4	2

Set Components

- Ureteral stent
- 2 fixation clamps
- Introducer made of polyurethane, green, approx. 45 cm long
- Flexible guidewire made of stainless steel with PTFE coating, fixed core, flexible safety tip and rigid tip: approx. 100 cm long in sets with cylindrical closed tip and approx. 150 cm long in sets with tip with central opening
- Sets for coaxial application include an additional ureteral catheter with central opening in stent size (see table)
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Maximum diameter of guidewire for paediatric Valve ureter stents with open tip

REF. 335403 / 335404 / 335405 / 335406 / 335407	
Ch. 3	Ch. 4
0.018"	0.028"



Valve DD ureteral stents

Made of polyurethane, with DD coupling and anti-reflux valve at the vesical end

Directable with anti-reflux valve

- Graduation in cm
- Radiopaque stent
- Pigtail and drainage eyes as indicated in the table
- Supplied as a preassembled unit
- Sterile, single use



Valve ureteral stents

Made of polyurethane, with anti-reflux valve at the vesical end

Anti-reflux valve


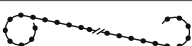
- Graduation in cm
- Radiopaque stent
- Pigtail and drainage eyes as indicated in the table
- Sterile, single use

When ordering, please always indicate Ref. and Order size

Integral ureteral stents

- Made of polyurethane
- Graduation in cm
  - Radiopaque stent
  - Pigtail and drainage eyes as indicated in the table
  - Sterile, single use

Integral ureteral stents

TIP AND EYE CONFIGURATION	REF. ACCORDING TO LENGTH OF STRAIGHT PART					ORDER SIZE (CH.)	QTY
	12 CM	15 CM	18 CM	20 CM	22 CM		
Closed tip - for simultaneous applications							
	334207	334231	334245	334232	334263	3/4/4.8	5
Closed tip - for coaxial applications							
	334700	334704	334705	334760	334706	3/4/4.8	5

Set components

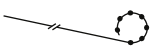
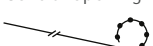
- Ureteral stent
- 2 fixation clamps
- Introducer made of polyurethane, green, approx. 45 cm long
- Flexible guidewire made of stainless steel with PTFE coating, fixed core, flexible safety tip and rigid tip: approx. 100 cm long in sets with cylindrical closed tip and approx. 150 cm long in sets with tip with central opening
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Maximum diameter of guidewire for paediatric integral ureteral stents with open tip

REF. 334700 / 334704 / 334705 / 334760 / 334706		
Ch. 3	Ch. 4	Ch. 4.8
0.018"	0.028"	0.028"

5.5 Special ureteral stents

Stents with one pigtail

TIP AND EYE CONFIGURATION	REF.	LENGTH OF STRAIGHT PART	ORDER SIZE (CH.)	QTY
Cylindrical closed 	334400	90 cm	4.8/6/7/8/9	6
Central opening 	334479	90 cm	4.8/6/7/8/9	6

Maximum diameter of guidewire for ureteral stents with one pigtail and open tip

REF. 334479				
Ch. 4.8	Ch. 6	Ch. 7	Ch. 8	Ch. 9
0.035"	0.035"	0.035"	0.038"	0.038"

Set components

- Ureteral stent
- Fixation clamp
- Flexible guidewire made of stainless steel with PTFE coating, fixed core, flexible safety tip and rigid tip, 150 cm long
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Ureteral stents with one pigtail

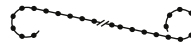
- Made of polyurethane
- Graduation in cm
  - Radiopaque stent
  - With one pigtail and drainage eyes only along the pigtail
  - Sterile, single use



Information

Each sales unit contains three stents with a marking ring at the proximal end and three stents without a marker. When stents have to be applied to both sides, the physician can place one stent with a marker and one without. The corresponding documentation in the patient card shows clearly which stent leads to which kidney..

Integral ureter stent for open surgery

TIP AND EYE CONFIGURATION	REF.	LENGTH OF STRAIGHT PART	ORDER SIZE (CH.)	QTY
	334209	26 cm	4.8/6/7	5
	334210	28 cm	4.8/6/7	

Set Components

- Ureteral stent
- Reinforced guidewire made of stainless steel with PTFE coating, 100 cm long
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging

Integral ureter stents for open surgery

- Made of polyurethane
- Graduation in cm
  - Radiopaque stent
  - With two cylindrical closed pigtails and drainage eyes as indicated in the table
  - Sterile, single use

When ordering, please always indicate Ref. and Order size

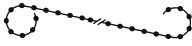
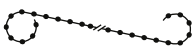


Integral ureteral stents

Made of polyurethane, short length

- Graduation in cm
- Radiopaque stent
- Pigtail and drainage eyes as indicated in the table
- Sterile, single use



Integral ureteral stent sets				
TIP AND EYE CONFIGURATION	REF.	LENGTH OF STRAIGHT PART	ORDER SIZE (CH.)	QTY
<div>Cylindrical closed</div> 	For simultaneous applications			
	334207	12 cm	7	5
<div>Central opening</div> 	For coaxial applications			
	334700	12 cm	6/7	5

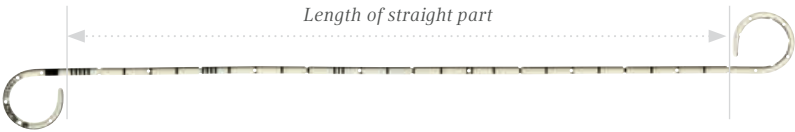
Maximum diameter of guidewire for short integral ureteral stents with open tip

REF. 334700	
Ch. 6	Ch. 7
0.035"	0.035"

Set components

- Ureteral stent
- 2 fixation clamps
- Introducer made of polyurethane, green, approx. 45 cm long
- Flexible guidewire made of stainless steel with PTFE coating, fixed core, flexible safety tip and rigid tip: approx. 100 cm long in sets with cylindrical closed tip and approx. 150 cm long in sets with tip with central opening
- Patient card
- 3 self-adhesive documentation labels
- Sterile packaging


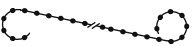

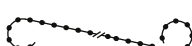
Length Measurement for Rüsch Ureteral Stents



5.6 ECO ureteral stent sets



ECO ureteral stents

TIP CONFIGURATION	EYE CONFIGURATION	REF.	LENGTH OF STRAIGHT PART	ORDER SIZE (CH.)	QTY
<div>Central opening</div> 		310000-026048 310000-026060 310000-026070 310000-026080	26 cm	4.8 6 7 8	10
		310000-028048 310000-028060 310000-028070 310000-028080	28 cm	4.8 6 7 8	10
		310000-030060 310000-030070	30 cm	6 7	10
		Set components <ul style="list-style-type: none"><li>• Ureteral stent</li><li>• <b>Without</b> guidewire</li><li>• Introducer made of polyurethane, transparent, approx. 45 cm long</li><li>• 3 self-adhesive documentation labels</li><li>• Sterile packaging</li></ul>			
<div>Cylindrical closed</div> 		310001-026048 310001-026060 310001-026070 310001-026080	26 cm	4.8 6 7 8	10
		310001-028048 310001-028060 310001-028070 310001-028080	28 cm	4.8 6 7 8	10
		Set components <ul style="list-style-type: none"><li>• Ureteral stent</li><li>• <b>With</b> guidewire made of stainless steel</li><li>• 2 fixation clamps</li><li>• Introducer made of polyurethane, transparent, approx. 45 cm long</li><li>• 3 self-adhesive documentation labels</li><li>• Sterile packaging</li></ul>			

ECO ureteral stent sets

Made of polyurethane

- Graduation in cm
- Radiopaque stent
- Pigtail and drainage eyes as indicated in the table
- Sterile, single use



Maximum diameter of guidewire for ECO ureteral stents with open tip


REF. 310000			
Ch. 4.8	Ch. 6	Ch. 7	Ch. 8
0.035"	0.035"	0.038"	0.038"

When ordering, please always indicate Ref. and Order size

# 6 Nephrostomy


## 6.1 Puncture sets

Superglide PCN 2-step puncture set							
	REF.	ORDER SIZE	PUNCTURE CANNULA		GUIDEWIRE O.D.	CATHETER EXTENDER I.D.	QTY
			O.D.	LENGTH			
	340012	Ch. 6	19.5 G / 0.95 mm	120 mm	0.028"	0.75 mm / 0.029"	1
		Ch. 8	17.5 G / 1.30 mm	200 mm	0.038"	1.00 mm / 0.039"	
		Ch. 10	17.5 G / 1.30 mm	200 mm	0.038"	1.00 mm / 0.039"	
	<b>Set components</b>						
	<ul style="list-style-type: none"><li>• 2-part puncture cannula with ultrasound marking at the tip (3 marking rings)</li><li>• Stiff guidewire with flexible J-curved tip in dispenser, with introducer, 800 mm long</li><li>• Catheter extender made of metal, 325 mm long, with plastic stylet</li><li>• Pigtail catheter made of polyurethane, radiopaque, approx. 30 cm long, central opening, hydrophilic coating, 6 drainage eyes, notches for retaining suture on the shaft, with Luer-Lock connector</li><li>• LL stopcock</li><li>• Adapter for urine collection bag (LL-funnel combination)</li><li>• 3 self-adhesive documentation labels</li><li>• Sterile, single use</li><li>• <b>Not made with natural rubber latex</b></li></ul>						

PCN 3-Step puncture set								
	REF.	ORDER SIZE	PUNCTURE CANNULA		GUIDEWIRE O.D.	DILATOR/SPLIT SHEATH		QTY
			O.D.	LENGTH		LENGTH	DILATOR O.D.	
	340013	Ch. 6	19.5 G / 0.95 mm	120 mm	0.028"	175 / 162.4 mm	Ch. 10	1
		Ch. 8	17.5 G / 1.30 mm	200 mm	0.038"	175 / 160.1 mm	Ch. 12	
		Ch. 10	17.5 G / 1.30 mm	200 mm	0.038"	175 / 157.8 mm	Ch. 14	
		<b>Set components</b> <ul style="list-style-type: none"><li>• 2-part puncture cannula with ultrasound marking at the tip (3 marking rings)</li><li>• Stiff guidewire with flexible J-curved tip in dispenser, with introducer, 800 mm long</li><li>• 2-part dilator, central opening, approx. 175 mm long, with split sheath, radiopaque</li><li>• Pigtail catheter made of polyurethane, radiopaque, approx. 30 cm long, central opening, 6 drainage eyes, notches for retaining suture on the shaft, with vertically slitted pigtail extender and Luer-Lock connector</li><li>• LL stopcock</li><li>• Adapter for urine collection bag (LL-funnel combination)</li><li>• 3 self-adhesive documentation labels</li><li>• Sterile, single use</li><li>• <b>Not made with natural rubber latex</b></li></ul>						

When ordering, please always indicate Ref. and Order size

## 6.2 Exchange sets

Superglide PCN exchange set					
	REF.	ORDER SIZE	GUIDEWIRE O.D.	CATHETER EXTENDER I.D.	QTY
	340014	Ch. 6	0.028"	0.75 mm / 0.029"	1
		Ch. 8	0.038"	1.00 mm / 0.039"	
		Ch. 10	0.038"	1.00 mm / 0.039"	
		<b>Set components</b> <ul style="list-style-type: none"><li>• Stiff guidewire with flexible J-curved tip in dispenser, with introducer, 800 mm long</li><li>• Catheter extender made of metal, 325 mm long, with plastic stylet</li><li>• Pigtail catheter made of polyurethane, radiopaque, approx. 30 cm long, central opening, hydrophilic coating, 6 drainage eyes, notches for retaining suture on the shaft, with Luer-Lock connector</li><li>• LL stopcock</li><li>• Adapter for urine collection bag (LL-funnel combination)</li><li>• 3 self-adhesive documentation labels</li><li>• Sterile, single use</li><li>• <b>Not made with natural rubber latex</b></li></ul>			



6.3 Catheters for percutaneous nephrostomy



Superglide balloon catheters

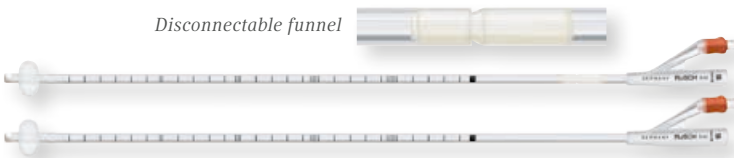
	REF.	ORDER SIZE	BALLOON	LENGTH	QTY
	Superglide nephrostomy balloon catheter, with hollow metal stylet				
	340408	Ch. 12/14	3ml	240mm	5
	340408	Ch. 16-24	5ml	240mm	
	340409	Ch. 12/14	3ml	410mm	
	340409	Ch. 16-24	5ml	410mm	

Superglide balloon catheters for nephrostomy


Made of silicone, with hollow metal stylet

- Tip with central opening, application over a guidewire of up to 0.038"
- Hydrophilic coating on first 5 cm
- 2 drainage openings in front of and 2 behind the balloon
- Valve for Luer and Luer-Lock syringe tips
- Transparent
- Cylindrical radiopaque tip
- Contrast stripe along the catheter shaft
- Graduation in cm
- 2 self-adhesive documentation labels
- Sterile, single use
- **Not made with natural rubber latex**

Disconnectable funnel



Balloon catheters

	REF.	ORDER SIZE	BALLOON	LENGTH	QTY
	With disconnectable funnel				
	340500	Ch. 12/14	3ml	450mm	5
		Ch. 16-24	5ml		
	With funnel which cannot be disconnected				
	340502	Ch. 12/14 Ch. 16-24	3ml 5ml	450mm	5

Balloon catheters for nephrostomy

Made of silicone, with the option of a funnel that is disconnectable

- Tip with central opening, radiopaque
- 2 drainage openings in front of and 2 behind the balloon
- Valve for Luer and Luer-Lock syringe tips
- Transparent
- Contrast stripe along the catheter shaft
- Graduation in cm
- Ref. 340502: exchange catheter for Ref. 340500
- 2 self-adhesive documentation labels
- Sterile, single use
- **Not made with natural rubber latex**

When ordering, please always indicate Ref. and Order size





7 Guidewires

Endoguide, stainless steel guidewires

- Smooth transition from the stiff wire to the flexible safety tip
- PTFE coating
- Rigid proximal end
- Distal tip straight or as J-curve
- In dispenser, with introducer
- Sterile, single use



Endoguide, stainless steel guidewires



	REF.	LENGTH	SIZE (O.D.)	QTY
	Endoguide stiff, straight tip			
	680141-000035	800mm	0.035"	5
	680141-000038		0.038"	
	680143-000028	1,000mm	0.028"	5
	680143-000035		0.035"	
	680143-000038		0.038"	
	680145-000028	1,200mm	0.028"	5
	680145-000035		0.035"	
	680145-000038		0.038"	
	Endoguide stiff, J-curved tip			
	680151-000028	800mm	0.028"	5
	680151-000035		0.035"	
	680151-000038		0.038"	
	680153-000028	1,000mm	0.028"	5
	680153-000035		0.035"	
	680153-000038		0.038"	
	680155-000028	1,200mm	0.028"	5
	680155-000035		0.035"	
	680155-000038		0.038"	

8 Stone extractors

8.1 Flat wire baskets

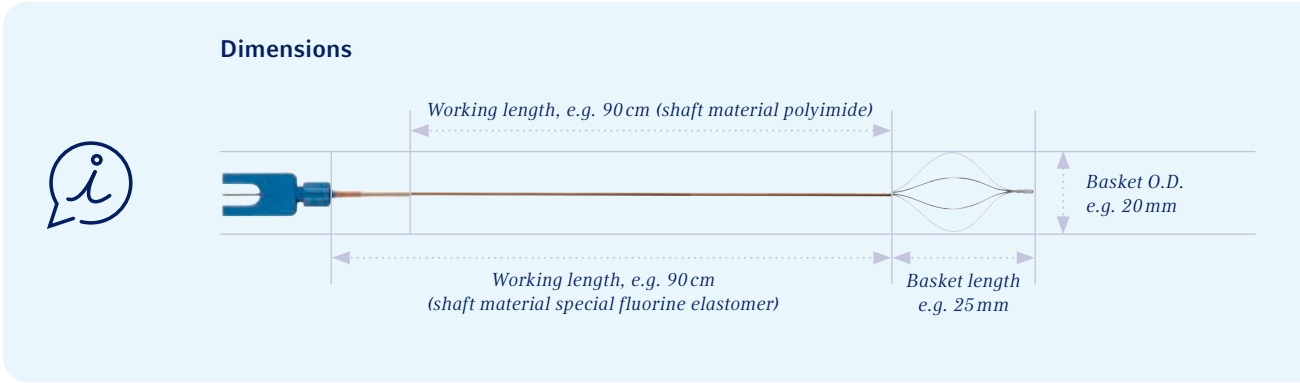


Rüsch-StoneX 4F

	REF. (QTY 1)				
	FOR TIP	ORDER SIZE Ø	WORKING LENGTH	BASKET O.D.	BASKET LENGTH
	<b>Rüsch-StoneX 4F with 4 flat wires and cylindrical tip</b>				
	342700	Ch. 2.5	90 cm	15 mm	20 mm
	342900	Ch. 2.5	115 cm	15 mm	20 mm
	342700	Ch. 3	90 cm	20 mm	25 mm
	342900	Ch. 3	115 cm	20 mm	25 mm
	<b>Rüsch-StoneX 4F with 4 flat wires and filiform tip made of polyurethane</b>				
	342800	Ch. 4	90 cm	20 mm	25 mm

Rüsch-StoneX flat wire baskets

- Form-stable flat wire basket made from stainless steel, symmetrical external shape
- Single-handed handgrip
- Clip-on retainer
- Radiopaque
- Packed in blister
- Shaft material: polyimide or special fluorine elastomer
- Sterile

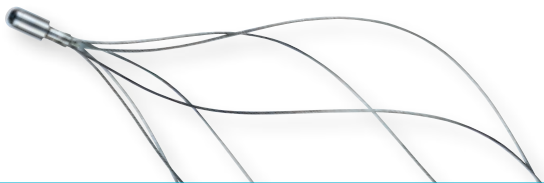


When ordering, please always indicate Ref. and Order size

8.2 Round wire baskets

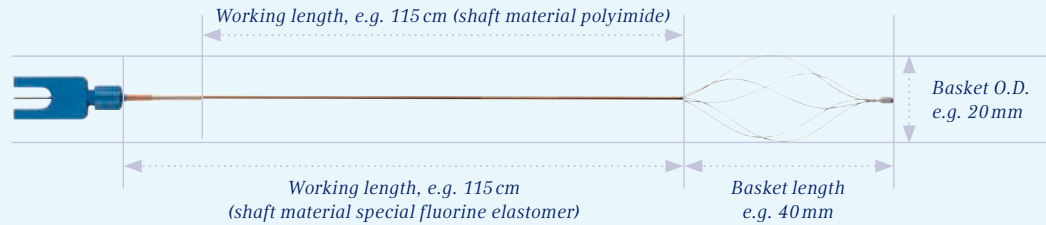
Rüsch-StoneX round wire baskets

- Form-stable round wire basket made from stainless steel, helical external shape
- Single-handed handgrip
- Clip-on retainer
- Shaft with cm graduations with 65 cm working length
- Radiopaque
- Packed in blister
- Shaft material: polyimide or special fluorine elastomer
- Sterile



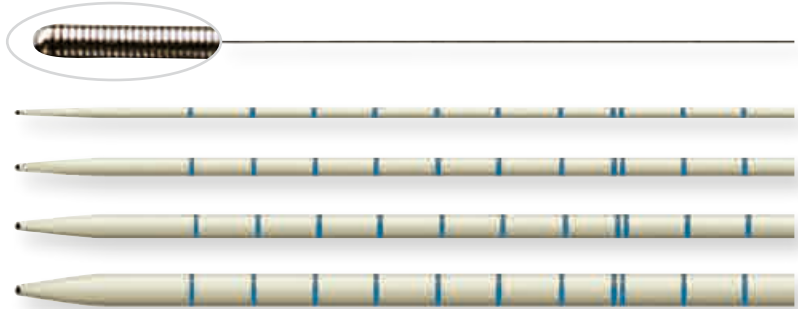
Rüsch-StoneX 4R, 5R, 6R					
	REF. (QTY 1)				
	FOR TIP	ORDER SIZE Ø	WORKING LENGTH	BASKET O.D.	BASKET LENGTH
	Rüsch-StoneX 4R with 4 round wires and cylindrical tip				
	341450	Ch. 2.5	90 cm	12 mm	20 mm
	341400	Ch. 2.5	115 cm	12 mm	20 mm
	341500	Ch. 3	65 cm	20 mm	40 mm
	341450	Ch. 3	90 cm	20 mm	40 mm
	341400	Ch. 3	115 cm	20 mm	40 mm
	341500	Ch. 4	65 cm	20 mm	40 mm
	341500	Ch. 5	65 cm	20 mm	40 mm
	Rüsch-StoneX 4R with 4 round wires and filiform tip made of polyurethane				
	341600	Ch. 4 Ch. 5	65 cm	20 mm	40 mm
	Rüsch-StoneX 5R with 5 round wires and cylindrical tip				
	341401	Ch. 3	115 cm	20 mm	40 mm
	341501	Ch. 4	65 cm		
	341501	Ch. 5	65 cm		
	Rüsch-StoneX 5R with 5 round wires and filiform tip made of polyurethane				
	341601	Ch. 4 Ch. 5	65 cm	20 mm	40 mm
	Rüsch-StoneX 6R with 6 round wires and cylindrical tip				
	341402	Ch. 3	115 cm	20 mm	40 mm
	341502	Ch. 4	65 cm		
	341502	Ch. 5	65 cm		
	Rüsch-StoneX 6R with 6 round wires and filiform tip made of polyurethane				
	341602	Ch. 4 Ch. 5	65 cm	20 mm	40 mm

Dimensions



9 Special products

9 Mechanical ureter dilators



Mechanical ureter dilators

Dilators made of polyurethane, guidewire made of stainless steel

- Polyurethane dilators
- Radiopaque
- Graduation in cm
- Sterile, single use

Mechanical ureter dilators

TIP CONFIGURATION	REF.	ORDER SIZE	LENGTH	QTY
 Conical with central opening	341000	Ureter dilator set, Vienna model, consisting of dilators in 9 different sizes and a guidewire with straight flexible safety tip and cm graduation up to 50 cm		1
 Flexible tip		Dilator set contains one size per each: Ch. 6/8/9/10/11/12/13/14/16	Approx. 60 cm	
		Guidewire Ø 0.037" (0.95 mm)	Approx. 150 cm	
 Conical with central opening	341001	Ureter dilator, individually packed		5
		Ch. 6/8/9/10/11/12/13/14/16	Approx. 60 cm	



# Index of reference numbers

REF.	PAGE
170003	13
170004	13
170605	4
170606	14
170630	7
171305	4
171505	4
172530	7
173230	9
173330	9
173420	14
173430	9
173800	10
173810	9
173830	9
173900	10
173910	9
173930	9
174000	10
174030	9
174100	10
174130	9
178000	4
178001	4
178003	13
178004	13
178305	4
180003	13
180105	5
180505	5
180605	5
180630	8
181305	5
181330	8
183230	12
183405	12

REF.	PAGE
183430	12
183530	8
183830	12
183930	12
183960	12
186005	5
186505	5
221800	16
223600	25
223602	25
223604	25
223606	25
252000	16
259700	16
263602	25
263702	25
263802	25
264001	26
264100	26
264102	26
264600	25
264701	24
264702	24
264710	24
264711	24
264810	24
264811	24
267500	27
310000	41
310001	41
331305	6
334200	35
334201	35
334202	35
334205	35
334206	35

REF.	PAGE
334207	38/40
334209	39
334210	39
334211	35
334212	35
334231	38
334232	38
334245	38
334247	30
334248	30
334249	30
334250	30
334263	38
334324	31
334326	31
334328	31
334400	39
334479	39
334700	38/40
334702	35
334703	35
334704	38
334705	38
334706	38
334710	35
334711	35
334714	35
334720	36
334721	36
334722	36
334723	36
334724	36
334740	30
334741	30
334742	30
334743	30

REF.	PAGE
334744	30
334745	30
334746	30
334751	31
334752	31
334753	31
334760	38
334764	36
334765	36
334766	36
334800	34
334801	34
334840	29
334841	29
334900	34
334901	34
334940	29
334941	29
335101	34
335140	29
335141	29
335142	29
335143	29
335166	34
335200	34
335201	34
335300	33
335301	33
335303	37
335304	37
335305	37
335306	37
335307	37
335309	37
335310	37
335311	33

REF.	PAGE
335312	37
335313	37
335314	37
335315	37
335316	37
335340	29
335341	29
335400	33
335401	33
335403	37
335404	37
335405	37
335406	37
335407	37
335411	33
335413	37
335414	37
335415	37
335416	37
335417	37
335441	29
335524	32
335526	32
335528	32
335620	32
335622	32
335624	32
335626	32
335628	32
340012	42
340013	42
340408	45
340409	45
340500	45
340502	45
340702	28

REF.	PAGE
340703	28
341000	49
341001	49
341400	48
341400	48
341401	48
341402	48
341450	48
341500	48
341501	48
341502	48
341600	48
341601	48
341602	48
342700	47
342800	47
342900	47
556000	15
660610	6
660630	7
661610	6
661630	7
662430	7
662530	7
662575	7
662630	7
662675	7
663230	11
663630	7
664130	11
664176	11
664230	11
664275	11
664330	11
664530	11
664575	11

REF.	PAGE
673044	19
673045	23
673049	19
673052	19
673053	22
673054	22
673064	21
673065	21
673081	20
673082	20
673083	20
680141	46
680143	46
680145	46
680151	46
680153	46
680155	46
850550	15
171100	18
335312	33

Teleflex is a global provider of medical technologies designed to improve the health and quality of people's lives. We apply purpose-driven innovation – a relentless pursuit of identifying unmet clinical needs – to benefit patients and healthcare providers. Our portfolio is diverse, with solutions in the fields of vascular and interventional access, surgical, anesthesia, cardiac care, urology, emergency medicine and respiratory care. Teleflex employees worldwide are united in the understanding that what we do every day makes a difference. For more information, please visit [teleflex.com](https://teleflex.com).

Teleflex is the home of Arrow, Deknatel, LMA, Pilling, QuikClot, Rüschi, UroLift, and Weck – trusted brands united by a common sense of purpose.

### Corporate Office

Phone +1 610 225 6800,  
550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

### Regional Offices

**United States:** Phone +1 919 544 8000, Toll Free 866 246 6990, [cs@teleflex.com](mailto:cs@teleflex.com), 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

**Latin America:** Phone +1 919 433 4999, [la.cs@teleflex.com](mailto:la.cs@teleflex.com), 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

**International:** Phone +353 (0)9 06 46 08 00, [orders.intl@teleflex.com](mailto:orders.intl@teleflex.com), Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

---

**Australia** 1300 360 226

**Austria** +43 (0)1 402 47 72

**Belgium** +32 (0)2 333 24 60

**Canada** +1 (0)800 387 9699

**China (Shanghai)** +86 (0)21 6060 7000

**China (Beijing)** +86 (0)10 6418 5699

**Czech Republic** +420 (0)495 759 111

**France** +33 (0)5 62 18 79 40

**Germany** +49 (0)7151 691 3004

**Greece** +30 210 67 77 717

**India** +91 (0)80 4093 4790

**Italy** +39 0362 58911

**Japan** +81 (0)3 6632 3600

**Korea** +82 2 536 7550

**Mexico** +52 55 5002 3500

**Netherlands** +31 (0)88 00 215 00

**New Zealand** 0800 601 100

**Poland** +48 22 4624032

**Portugal** +351 225 491 051

**Singapore** (SEA non-direct sales countries) +65 6439 3000

**Slovak Republic** +421 (0)3377 254 28

**South Africa** +27 (0)11 807 4887

**Spain** +34 91 198 84 31

**Switzerland** +41 (0)31 818 40 90

**United Kingdom** +44 (0)1494 53 27 61

**For more information, please visit [teleflex.com](https://teleflex.com)**

---

Teleflex, the Teleflex logo, Rüschi, AquaFlate, Endoguide, Rüschi Brilliant, Rüschi Gold, Rüschi StoneX, SilFlate, Silkospray, SoftSimplastich, Superglide and Wiruthan are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. © 2022 Teleflex Incorporated. All rights reserved. Revised: 01/2022.

MCI-100224 · REV 0 · 01 22 PDF

**Teleflex**